

Mouse Myeloperoxidase ELISA Instructions

CAT. No: AEM0030

CONTENT

No.	Content	Cat. No	Volume
1	CP (Cotaed Plate)	EM0030 CP	96 wells
2	S (Standard)	EM0030S,S1~ S7,S0	9 vials
3	DA (Detect Antibody)	EM0030 DA	6 ml/bottle
4	SH(Streptavidin-HRP)	ESH03	12 ml/bottle
5	AB(Assay Buffer 1×)	EAB01	12 ml/bottle
6	TS (TMB Substrate)	ETS01	12 ml/bottle
7	SS (Stop Solution)	ESS01	12 ml/bottle
8	WB (Wash Buffer 10×)	EWB01	50 ml/bottle
9	SF (Sealer Film)	ESF01	6 pieces

NOTE: After the kit is opened, the stabilization period of each content is 30 days, so please use the kit within 30 days after opening.


REAGENT PREPARATION

Washing Buffer (1×) Preparation: Pour entire contents (50 ml) of the **Washing Buffer Concentrate** (10×) into a clean 500 ml graduated cylinder. Bring to final volume of 500 ml with glass-distilled or deionized water. Transfer to a clean wash bottle and store at 2 to 25°C.

Standard Curve Preparation: S1 to S7 and S0 is ready to use for serum and plasma.

Other sample type, prepare the standard curve with whatever buffer (SPB, Sample Prepared Buffer) is used to prepare the sample, such as cell culture supernatant, tissue grinding liquid, cell lysate, etc. Urine sample use AB (Assay Buffer) prepare standard curve.

The Mouse Myeloperoxidase Standard EM0030S 1000 ng/ml 30 µl + 270 µl SPB serves as the high standard (100 ng/ml). Pipette 200 µl of SPB into each tube. Use the high standard to produce a 1:2 dilution series. Mix each tube thoroughly before the next transfer. SPB serves as the zero standard (0 pg/ml).

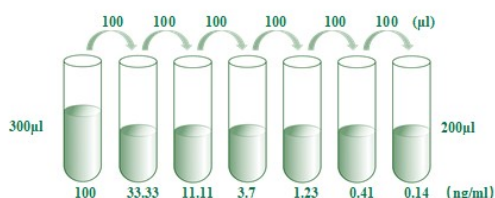


ASSAY PROCEDURE

Bring all reagents and samples to room temperature before use.

1. Prepare all reagents and working standards as directed in the previous sections.
2. Remove excess **CP** (Coated Plate) strips from the plate frame, return them to the foil pouch and reseal.
3. Add 50 µl of **AB** (Assay Buffer) to each well.
4. Add 10 µl of **Standard or sample** per well. Ensure reagent addition is uninterrupted and completed within 15 minutes.
5. Add 50 µl of **DA** (Detect Antibody) to each well.
6. Cover with an **SF** (Sealer Film). Incubate at room temperature (18 to 25°C) for 1 hour on a microplate **shaker** set at 500 rpm.
7. Aspirate each well and **wash**, repeating the process four times. Wash by filling each well with WB (Washing Buffer 300 µl). Complete removal of liquid at each step is essential to good performance. After the last wash, remove any remaining WB (Washing Buffer) by aspirating or decanting. Invert the plate and **blot** it against clean paper towels.
8. Add 100 µl of **SH** (Streptavidin-HRP) to each well.
9. Cover with a new **SF** (Sealer Film). Incubate at room temperature (18 to 25°C) for 30 minutes on a microplate **shaker** set at 500 rpm.
10. Repeat aspiration/**wash** as in step 7.
11. Add 100 µl of **TS** (TMB Substrate) to each well. Incubate for 5-30 minutes at room temperature.
12. Add 100 µl of **SS** (Stop Solution) to each well.
13. Determine the optical density within 30 minutes, using microplate **reader** set to 450 nm corrected with 570 nm or 630 nm.

Mouse Myeloperoxidase Typical Standard

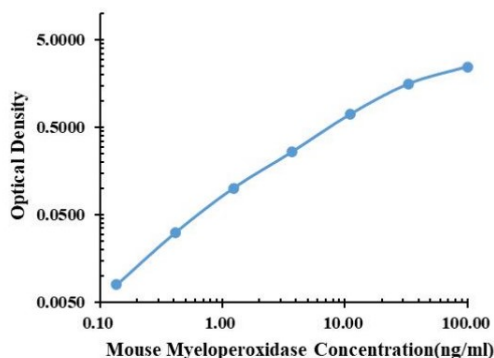

ng/ml	O.D.		Average	Corrected
0.00	0.0213	0.0192	0.0203	
0.14	0.0285	0.0280	0.0283	0.0080
0.41	0.0525	0.0505	0.0515	0.0312
1.23	0.1142	0.1285	0.1214	0.1011
3.70	0.2889	0.2755	0.2822	0.2619
11.11	0.7461	0.7085	0.7273	0.7070
33.33	1.6370	1.5500	1.5935	1.5732
100.00	2.5460	2.4060	2.4760	2.4557

SENSITIVITY

The minimum detectable dose (MDD) of mouse Myeloperoxidase is typically less than 0.10 pg/ml.

The MDD was determined by adding two standard deviations to the mean optical density value of ten zero standard replicates and calculating the corresponding concentration.

PRECISION

- Intra-assay Precision (Precision within an assay) Three samples of known concentration were tested twenty times on one plate to assess intra-assay precision.
- Inter-assay Precision (Precision between assays)

	Intra-assay Precision				Inter-assay Precision		
	S1	S2	S3		S1	S2	S3
Sample Number	22	22	22		6	6	6
Average (pg/ml)	9.6	3.7	1.1		10.0	4.1	1.1
Standard deviation	0.7	0.3	0.1		0.6	0.2	0.1
Coefficient of variation (%)	6.9	7.0	8.0		6.5	5.8	5.8

RECOVERY

The spike recovery was evaluated by spiking 3 levels of mouse Myeloperoxidase into health mouse serum sample. The un-spiked serum was used as blank in this experiment.

The recovery ranged from 100% to 129% with an overall mean recovery of 117% .

LINEARITY

To assess the linearity of the assay, five samples were spiked with high concentration of Myeloperoxidase in mouse serum and diluted with Sample Diluent to produce samples with values within the dynamic range of the assay.

The linearity ranged from 101% to 127% with an overall mean recovery of 115%.

SAMPLE VALUES

Serum/Plasma - Thirty samples from apparently healthy mice were evaluated for the presence of Myeloperoxidase in this assay. No medical histories were available for the donors.

Sample Matrix	Sample Evaluated	Range (pg/ml)	Detectable %	Mean of Detectable (pg/ml)
Serum	30	31.55-52.87	100	45.23

n.d. = non-detectable. Samples measured below the sensitivity are considered to be non-detectable.

DISCLAIMER AND VERSION

For research use only. Not for use in clinical diagnostic procedures.

Copyright 2024 AREX Biosciences Ltd. Version.240820